



Reciprocal Shaker

Reciprocal Shaker

The reciprocal shaker features a left to right linear action to provide gentle to vigorous shaking for a variety of applications. The operation may be continuous or timed, with the integral electronic timer ensuring accurate repeatability of time-sensitive incubations. Speed and time setting are straightforward in the LED display. Its load weight is up to 15kg. An audible signal accompanies automatic switch-off to indicate the completion of a set time period. The interchangeable / stacking platforms, and accessories enable a wide choice of vessels to be used, including bottles, flasks and beakers, dishes, boxes and petri-dishes.


Features

- >Continuous or timed operation with automatic switch-off
- >Variable shaking speed from 5 to 100rpm
- >Interchangeable / stacking platforms, and accessories for a variety of vessels
- >15kg carry capability
- >26 x 26cm platform & 30 x 30cm platform available



CE

Cat. No. :MS-NRC-30

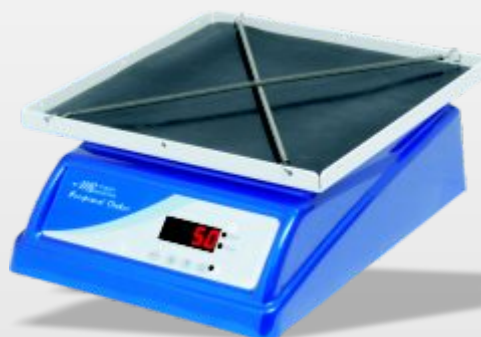

CE

Cat. No. :MS-NRC




Specification	MS-NRC	MS-NRC-30
Motion	Reciprocal	
Maximum Stroke Length	20mm	
Controller	Digital microprocessor controller	
Speed / Inc	5 ~ 100rpm / 1rpm	
Timer	1~9999 mins with alarm, continuous / 1 min	
Carry Capability	15kg	
Motor	DC brushless type	
Operation Temperature	Ambient to 40 °C	
Platform Dimension (W x L)	26 x 26cm	30 x 30cm
Platform Material	Painted aluminum alloy	Painted iron metal
Stacking Platform	Yes	
Unit Dimension (W x L x H)	260 x 310 x 130mm	
Construction	Frame retard PC and painted iron metal rare and bottom plates	
Weight	approx. 7.0kg	
Rated Voltage	110V / 220V selectable	



Ordering Information

MS-NRC	MS Reciprocal Shaker with 26 x 26cm platform and flat non-slip rubber mat
MS-NRC-FH250	MS-NRC with 5 of 250ml flask holder
MS-NRC-FH500	MS-NRC with 4 of 500ml flask holder
MS-NRC-30	MS-NRC with 30 x 30cm platform and flat non-slip rubber mat
MS-NRC-30-FH250	MS-NRC-30 with 9 of 250ml flask holder
MS-NRC-30-FH500	MS-NRC-30 with 9 of 500ml flask holder
MS-SP	Strip springs (pack of 2)
MS-P2626	Additional 26 x 26cm platform with 4 of adjustable 10cm or 20cm pillars and flat non-slip rubber mat
MS-P3030	Additional 30 x 30cm platform with 4 of adjustable 10cm or 20cm pillars and flat non-slip rubber mat
MS-DIMPLED	26 x 26cm dimpled mat
MS-DIMPLED-30	30 x 30cm dimpled mat

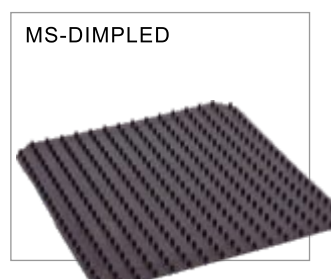



Contact Information

www.majorsci.com

info@majorsci.com

19959 Sea Gull Way
Saratoga, CA 95070
U.S.A

T/1-408-366-9866

F/1-408-446-1107

Innovative Life Sciences Tools