



50
YEARS



 **Texwipe®**
An ITW Company

Cleanroom Swabs

Microdenier | Polyester | Foam | ESD | Cotton | Specialty

Sterile | Sterile Dry Collection and Transport System | Cleaning Validation & TOC Validation Kits
Alpha® | Absorbond® | CleanFoam® | ESD-Safe | Cotton | Specialty | General Purpose

Table of Contents

	SterilePg. 2
	Sterile Dry Collection and Transport System. .Pg. 3
	Cleaning Validation . . .Pg. 4
	Cleaning Validation Kits Pg. 4
	MicrodenierPg. 5
	PolyesterPg. 6
	CleanFoam®.Pg. 7
	ESD-SafePg. 8
	Cotton.Pg. 9
	SpecialtyPg. 10
	General Purpose . . . Pg. 11

Our highly professional team will help you select and develop products for any critical environment application. Contact information is located on the back of this brochure or visit us at www.texwipe.com.

About Texwipe

Texwipe's commitment to innovation, leadership and quality in cleanroom consumable products spans more than fifty years. We invest in technology to respond to our customers' evolving needs in contamination control.



Texwipe's Core Values are:

- **Innovation** – Texwipe pioneers the latest technologies to provide innovation in contamination control products and processes.
- **Quality** – Texwipe product quality is maintained by the most advanced testing and quality control standards in the industry.
- **Technology Leadership** – Texwipe leads our industry in testing metrology, methods and processes to reduce contamination.

Throughout Texwipe's global operations, we support our customers with products designed to exceed the requirements for cleanroom consumable products. Our highly professional team will help you select and develop products for any critical environment application.

Texwipe CleanTips®

Texwipe Cleanroom Swabs, with their trademark green handles, are the standard in cleanrooms across the globe. Engineered for superior cleanliness and consistency for a variety of applications, Texwipe has a swab that will provide the best contamination performance available.

Our products range from Microdenier, Polyester, CleanFoam®, ESD-Safe, General Purpose, Cotton and our Specialty line featuring the pre-wetted CrushTube™ Swab.

Texwipe Cleanroom Swabs are manufactured with a complete thermal bond, eliminating adhesive contamination*. Our swabs are manufactured using high-precision automated processes, resulting in swabs constructed to exacting and consistent tolerances, lot coded for traceability and quality control.

Custom Swabs

At Texwipe® we understand that the best swab for your application may be something unique. We can create a custom swab for your specific applications.

Sterile Swabs

We can create sterile versions of all Texwipe swabs or something unique and specific to your application.

*excludes Cotton Swabs



www.texwipe.com

Sterile Swabs

Texwipe's Sterile Swabs are available in USP-grade cotton, polyurethane foam and medical-grade spun polyester. The swabs are made to exacting and consistent tolerances using high-precision automated processes. The easy-to-open, individually wrapped swabs ensure sterility at point-of-use. The swabs are visible through clear packaging which also offers visibility of seal integrity.



Cotton Series

Part No. STX705P	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	7.0mm (0.276")	18.2mm (0.717")	2.5mm (0.098")	2.5mm (0.098")	134.0 mm (5.276")	155.0mm (6.102")	1 swab per sleeve 50 swabs per box 10 boxes per case




Part No. STX705W	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	7.0mm (0.276")	18.2mm (0.717")	2.5mm (0.098")	2.5mm (0.098")	134.0mm (5.276")	155.0mm (6.102")	1 swab per sleeve 50 swabs per box 10 boxes per case




Foam Series

Part No. STX708A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	19.7mm (0.776")	9.4mm (0.370")	26.4mm (1.040")	6.6mm (0.260")	3.2mm (0.126")	105.8mm (4.165")	132.2mm (5.205")	1 swab per sleeve 50 swabs per box 10 boxes per case

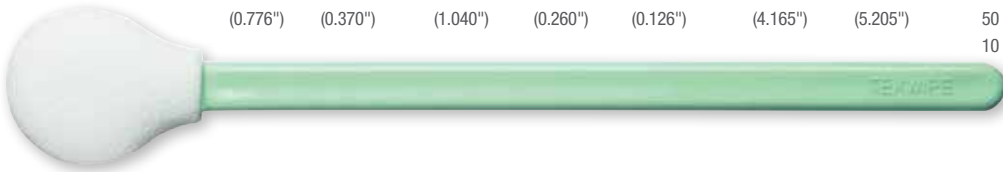

Part No. STX712A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	12.7mm (0.500")	7.8mm (0.307")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	1 swab per sleeve 50 swabs per box 10 boxes per case

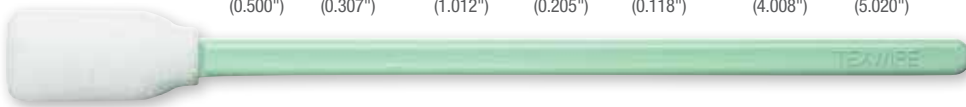

Polyester Series

Part No. STX763	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	5.8mm (0.228")	5.8mm (0.228")	18.0mm (0.709")	2.5mm (0.098")	2.5mm (0.098")	134.0mm (5.276")	155.0mm (6.102")	1 swab per sleeve 50 swabs per box 10 boxes per case



Part No. STX764	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	4.6mm (0.181")	4.6mm (0.181")	17.0mm (0.669")	2.5mm (0.098")	2.5mm (0.098")	134.0mm (5.276")	155.0mm (6.102")	1 swab per sleeve 50 swabs per box 10 boxes per case



Sterile Dry Collection and Transport System

Texwipe's Sterile Dry Collection and Transport System Swabs consist of sterile ready-to-use units for collection and transport of bacteriological specimens. The system is designed to maintain integrity of both sample collection and sample transport to the lab. Swabs are constructed from 100% USP-grade spun cotton or 100% USP-grade spun polyester. The swabs are bonded securely to a polystyrene shaft by an aqueous-based adhesive. 6.1" (155 mm) long swabs enable access to hard-to-reach sampling sites. The swabs are made to exacting and consistent tolerances using high-precision automated processes.



Cotton Series

Part No. STX705PT	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	7.0mm (0.276")	18.2mm (0.717")	2.5mm (0.098")	2.5mm (0.098")	134.0 mm (5.276")	155.0mm (6.102")	1 swab per tube 100 tubes per box 5 boxes per case



Polyester Series

Part No. STX764T	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	4.6mm (0.181")	4.6mm (0.181")	17.0mm (0.669")	2.5mm (0.098")	2.5mm (0.098")	134.0 mm (5.276")	155.0mm (6.102")	1 swab per tube 100 swabs per box 5 boxes per case



Sterile Available
All Texwipe swabs available in sterile upon request.

Texwipe is certified to ISO 9001 and ISO 13485.


Cleaning Validation Swabs

Engineered to deliver consistent and accurate sample recovery for cleaning validation for both specific (HPLC, UV-Vis, IMS) and nonspecific (TOC) analytical methods. The swabs are effective with a wide range of diluents.


Alpha®

Alpha series swab additionally processed to <50 µg/L TOC per swab.


Part No. TX714K	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	20 swabs per bag 50 bags per case

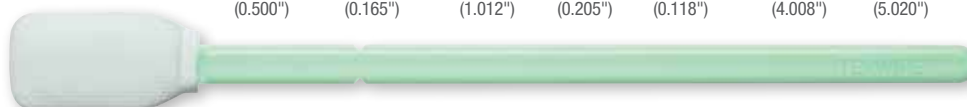

Part No. TX715	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case

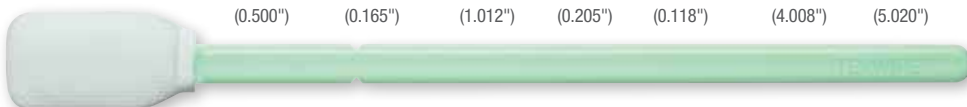


Part No. TX761K	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	6.8mm (0.268")	2.8mm (0.110")	16.8mm (0.661")	3.2mm (0.126")	3.2mm (0.126")	145.5mm (5.728")	162.3mm (6.390")	20 swabs per bag 50 bags per case

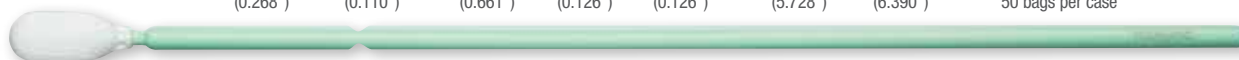

Absorbond®

Surfactant-free Absorbond Series swab for HPLC/IMS analysis.

Part No. TX716	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	20 swabs per bag 50 bags per case

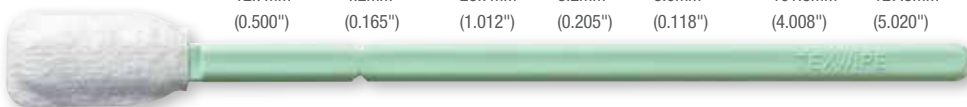

TOC Cleaning Validation Kits

Texwipe's TOC Cleaning Validation Kits are designed to simplify sampling as part of a cleaning validation protocol. The swabs and vials are further cleaned using a proprietary processes to make these kits excellent for use in Total Organic Carbon (TOC) analysis.

The kits are an organized, convenient and secure method for transporting swabs and vials from storage area to sampling sites to laboratory with minimal chance of contamination. Packaged in cleanroom-compatible, recyclable polypropylene boxes.

Available with Texwipe Cleanroom Swab TX714K (TX3340 and TX3342) or TX761K (TX3343).



TX3340
TOC Cleaning
Validation Kit

Product Number	Description	Packaging
TX3340	TOC Cleaning Validation Kit: 12 vials + 24 swabs (TX714K) + 12 labels	18 kits/case
TX3342	TOC Cleaning Validation Kit, Bulk: 72 vials + 144 swabs (TX714K) + 72 labels	1 kit/case
TX3343	TOC Cleaning Validation Kit: 12 vials + 24 swabs (TX761K) + 12 labels	18 kits/case

Polyester Swabs

Microdenier Series

Texwipe's Microdenier Swab Series is made from the cleanest microdenier material, 100% polyester, nylon-free. The complete thermal bond construction eliminates adhesive contamination.

Cleanroom laundered, providing ultra low levels of particles, NVRs (non volatile residues) and ions, and engineered to provide enhanced absorbency and quick wicking performance. Streak-free cleaning with scratch-sensitive surfaces and compatible with QAC (quaternary ammonium chloride) cleaners and disinfectants. 100% virgin polypropylene handle.


Cleanroom manufactured, made to exacting and consistent tolerances using high-precision automated processes. Lot coded for traceability and quality control.

Packaged in a reclosable, silicone-free bag.



Reclosable, silicone-free bag

Part No. TX714MD	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	100 swabs per bag 8 bags per case

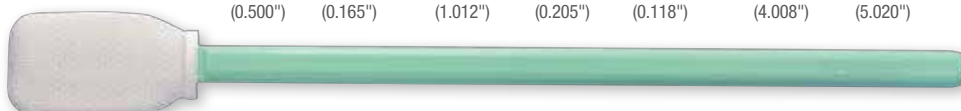

Part No. TX758MD	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.2mm (0.126")	2.5mm (0.098")	11.0mm (0.433")	2.2mm (0.087")	2.2mm (0.087")	60.0mm (2.362")	70.0mm (2.756")	500 swabs per bag 5 bags per case



Part No. TX761MD	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	6.8mm (0.268")	3.6mm (0.142")	16.8mm (0.661")	3.2mm (0.126")	3.2mm (0.126")	145.5mm (5.728")	162.5mm (6.398")	100 swabs per bag 8 bags per case




Sterile Available
All Texwipe swabs available in sterile upon request.

Polyester Swabs

Alpha® Series


Texwipe's Alpha® Swab Series is made from the cleanest 100% polyester knit materials with a double layer of fabric for enhanced absorbency. A complete thermal bond construction eliminates adhesive contamination. Cleanroom laundered and manufactured, providing ultra low levels of contaminants. Engineered to resist snagging and abrading. 100% virgin polypropylene handle. Made to exacting, consistent tolerances using high-precision automated processes. Packaged in a silicone-free, amide-free bag.

Part No. TX714A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
------------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

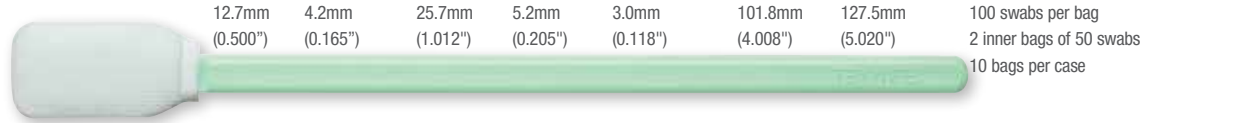


12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case
--------------------	-------------------	--------------------	-------------------	-------------------	---------------------	---------------------	---

Part No. TX743B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
------------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

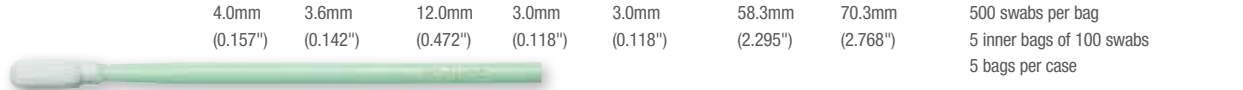


4.0mm (0.157")	3.6mm (0.142")	12.0mm (0.472")	3.0mm (0.118")	3.0mm (0.118")	58.3mm (2.295")	70.3mm (2.768")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case
-------------------	-------------------	--------------------	-------------------	-------------------	--------------------	--------------------	---

Part No. TX754B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
------------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

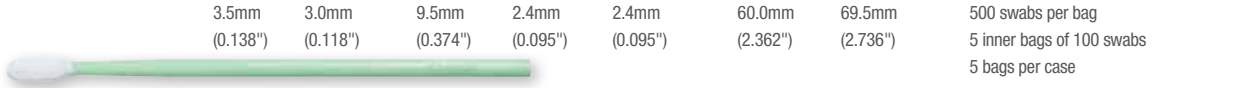


3.5mm (0.138")	3.0mm (0.118")	9.5mm (0.374")	2.4mm (0.095")	2.4mm (0.095")	60.0mm (2.362")	69.5mm (2.736")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case
-------------------	-------------------	-------------------	-------------------	-------------------	--------------------	--------------------	---

Part No. TX758B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
------------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

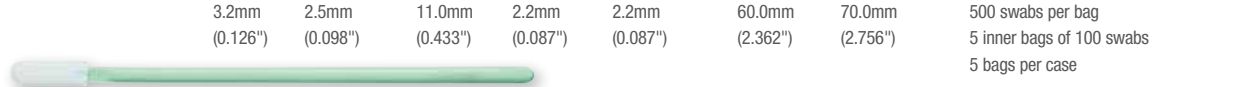

3.2mm (0.126")	2.5mm (0.098")	11.0mm (0.433")	2.2mm (0.087")	2.2mm (0.087")	60.0mm (2.362")	70.0mm (2.756")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case
-------------------	-------------------	--------------------	-------------------	-------------------	--------------------	--------------------	---

Part No. TX761	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
-----------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

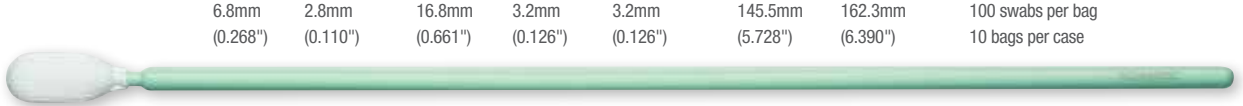


6.8mm (0.268")	2.8mm (0.110")	16.8mm (0.661")	3.2mm (0.126")	3.2mm (0.126")	145.5mm (5.728")	162.3mm (6.390")	100 swabs per bag 10 bags per case
-------------------	-------------------	--------------------	-------------------	-------------------	---------------------	---------------------	---------------------------------------

Absorbond® Series

Texwipe's Absorbond® Swabs Series is made from 100% polyester (hydroentangled) nonwoven material double layer for enhanced absorbency. The complete thermal bond construction eliminates adhesive contamination. Cleanroom processed and manufactured, providing ultra low levels of contamination control. Engineered for use on sensitive surfaces. 100% virgin polypropylene handle. Made to exacting, consistent tolerances using high-precision automated processes. Packaged in a silicone-free, amide-free bag.

Part No. TX716	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
-----------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

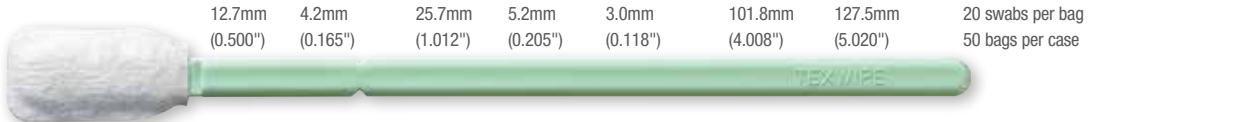


12.7mm (0.500")	4.2mm (0.165")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	20 swabs per bag 50 bags per case
--------------------	-------------------	--------------------	-------------------	-------------------	---------------------	---------------------	--------------------------------------

Part No. TX759B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
------------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------



3.2mm (0.126")	1.5mm (0.059")	10.0mm (0.394")	2.2mm (0.087")	2.2mm (0.087")	60.0mm (2.362")	70.0mm (2.756")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case
-------------------	-------------------	--------------------	-------------------	-------------------	--------------------	--------------------	---

Part No. TX762	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
-----------------------	------------	----------------	-------------	--------------	------------------	---------------	-------------------	-----------

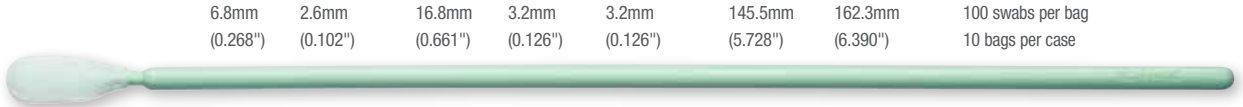

6.8mm (0.268")	2.6mm (0.102")	16.8mm (0.661")	3.2mm (0.126")	3.2mm (0.126")	145.5mm (5.728")	162.3mm (6.390")	100 swabs per bag 10 bags per case
-------------------	-------------------	--------------------	-------------------	-------------------	---------------------	---------------------	---------------------------------------

CleanFoam[®] Swabs

CleanFoam[®] Series A

Texwipe's CleanFoam[®] Swab Series is made with 100% polyurethane foam, 100 pores per inch, for good chemical compatibility with a variety of solutions. The complete thermal bond construction eliminates adhesive contamination.

100% virgin polypropylene handle. Cleanroom processed and made to exacting, consistent tolerances using high-precision automated processes. Packaged in a silicone-free, amide-free bag.


Foam Series A features 100 swabs per bag.




Part No. TX706A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	10.9mm (0.429")	8.5mm (0.335")	27.0mm (1.063")	6.3mm (0.248")	6.3mm (0.248")	54.0mm (2.126")	81.0mm (3.189")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case



Part No. TX707A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	15.0mm (0.591")	8.8mm (0.346")	25.0mm (0.984")	6.6mm (0.260")	3.2mm (0.126")	106.0mm (4.173")	131.0mm (5.157")	100 swabs per bag; 2 inner bags of 50 swabs; 10 bags per case

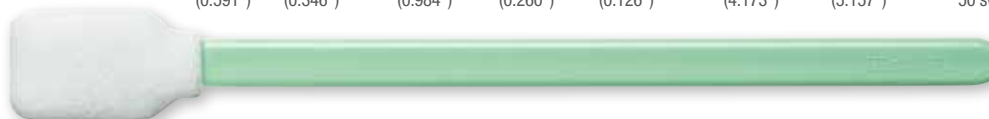

Part No. TX708A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
<i>Sterile STX708A</i>	19.7mm (0.776")	9.4mm (0.370")	26.4mm (1.039")	6.6mm (0.260")	3.2mm (0.126")	105.8mm (4.165")	132.2mm (5.205")	100 swabs per bag; 2 inner bags of 50 swabs; 10 bags per case

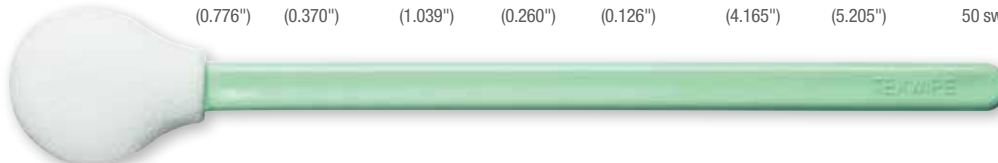

Part No. TX709A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	4.8mm (0.189")	4.2mm (0.165")	12.0mm (0.472")	4.3mm (0.169")	4.3mm (0.169")	92.5mm (3.642")	104.0mm (4.094")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case



Part No. TX710A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	5.9mm (0.232")	20.4mm (0.803")	4.3mm (0.169")	4.3mm (0.169")	92.5mm (3.642")	113.0mm (4.449")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case



Part No. TX712A	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
<i>Sterile STX712A</i>	12.7mm (0.500")	7.8mm (0.307")	25.7mm (1.012")	5.2mm (0.205")	3.0mm (0.118")	101.8mm (4.008")	127.5mm (5.020")	100 swabs per bag 2 inner bags of 50 swabs 10 bags per case



Sterile Available

All Texwipe swabs available in sterile upon request.

Texwipe is certified to ISO 9001 and ISO 13485.

CleanFoam® Swabs

CleanFoam® Series B

Texwipe's CleanFoam® Swab Series is made with 100% polyurethane foam, 100 pores per inch, for good chemical compatibility with a variety of solutions. The complete thermal bond construction eliminates adhesive contamination.

100% virgin polypropylene handle. Cleanroom processed and made to exacting, consistent tolerances using high-precision automated processes. Packaged in a silicone-free, amide-free bag.

Foam Series B features 500 swabs per bag.



Part No. TX740B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	6.2mm (0.244")	5.0mm (0.197")	17.0mm (0.669")	3.2mm (0.126")	3.2mm (0.126")	146.0mm (5.748")	163.0mm (6.417")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX741B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.4mm (0.134")	3.0mm (0.118")	10.5mm (0.413")	2.4mm (0.094")	2.4mm (0.094")	56.8mm (2.236")	67.3mm (2.650")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX742B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.6mm (0.142")	3.5mm (0.138")	12.0mm (0.472")	3.0mm (0.118")	3.0mm (0.118")	59.0mm (2.323")	71.0mm (2.795")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX751B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.5mm (0.138")	3.5mm (0.138")	16.0mm (0.630")	2.4mm (0.094")	2.4mm (0.094")	51.0mm (2.008")	67.0mm (2.638")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX752B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.4mm (0.134")	3.4mm (0.134")	20.0mm (0.787")	2.5mm (0.098")	2.5mm (0.098")	50.0mm (1.969")	70.0mm (2.756")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX757B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.2mm (0.126")	2.6mm (0.102")	10.0mm (0.395")	2.2mm (0.087")	2.2mm (0.087")	59.7mm (2.350")	69.7mm (2.744")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



ESD-Safe Swabs

Texwipe's ESD Swab Series features a Stat-Rite® handle, a patented Inherently Dissipative Polymer, to provide clean ESD protection without contributing particulate contaminants to the process or environment. The complete thermal bond construction eliminates adhesive contamination. Packaged in an anti-static bag.

Available in a choice of CleanFoam®, Alpha® Polyester, or Stick.



CleanFoam® ESD


Part No. TX750E	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.5mm (0.138")	3.5mm (0.138")	10.0mm (0.394")	3.0mm (0.118")	3.0mm (0.118")	78.5mm (3.090")	88.5mm (3.484")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX753E	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.7mm (0.146")	3.4mm (0.134")	10.0mm (0.394")	3.0mm (0.118")	3.0mm (0.118")	80.0mm (3.150")	90.0mm (3.543")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX757E	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.8mm (0.150")	3.4mm (0.134")	10.0mm (0.394")	3.0mm (0.118")	3.0mm (0.118")	80.7mm (3.177")	90.7mm (3.571")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case

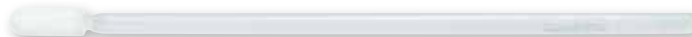

Alpha® Polyester ESD

Part No. TX758E	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	4.0mm (0.157")	3.0mm (0.118")	10.0mm (0.394")	3.0mm (0.118")	3.0mm (0.118")	80.5mm (3.169")	90.5mm (3.563")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Part No. TX761D	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	6.8mm (0.268")	4.0mm (0.157")	17.0mm (0.669")	3.0mm (0.118")	3.0mm (0.118")	147.0mm (5.787")	164.0mm (6.457")	100 swabs per bag 10 bags per case



Stick ESD

Part No. TX769E	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	1.0mm (0.039")	1.0mm (0.039")	9.6mm (0.378")	2.4mm (0.094")	2.4mm (0.094")	56.5mm (2.224")	66.1mm (2.602")	500 swabs per bag 5 inner bags of 100 swabs 5 bags per case



Sterile Available

All Texwipe swabs available in sterile upon request.

Texwipe is certified to ISO 9001 and ISO 13485.

Cotton Swabs

Texwipe's Cotton and Cotton/Foam Swab Series is designed to maximize the characteristics of cotton. Available in 100% cotton or cotton/polyurethane foam. These swabs are bonded securely to a wooden shaft by an aqueous-based adhesive.

Part No. TX700B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	13.5mm (0.531")	13.5mm (0.531")	28.0mm (1.102")	2.5mm (0.098")	2.5mm (0.098")	126.0mm (4.961")	154.0mm (6.063")	500 swabs per bag 10 inner bags of 50 swabs 5 bags per case




Part No. TX705	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	7.0mm (0.276")	17.0mm (0.669")	2.5mm (0.098")	2.5mm (0.098")	135.0mm (5.315")	152.0mm (5.984")	1,000 swabs per bag 10 inner bags of 100 swabs 5 bags per case



Part No. TX705P	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	7.0mm (0.276")	7.0mm (0.276")	17.0mm (0.669")	2.5mm (0.098")	2.5mm (0.098")	135.0mm (5.315")	152.0mm (5.984")	1,000 swabs per bag 10 inner bags of 100 swabs 5 bags per case



Part No. TX720B	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	11.8mm (0.465")	10.0mm (0.394")	26.3mm (1.035")	2.5mm (0.098")	2.5mm (0.098")	127.0mm (5.000")	154.0mm (6.063")	500 swabs per bag 10 inner bags of 50 swabs 5 bags per case

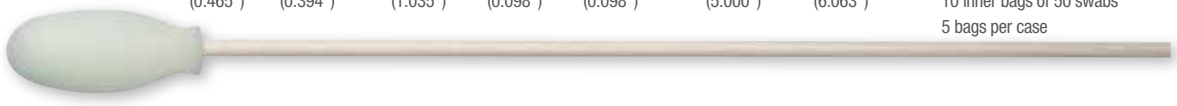

Specialty Swabs & Tools

TX726 is our innovative CrushTube Swab, featuring 91% IPA in an enclosed tube. Packaged in a resealable bag.

Part No. TX726	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	8.3mm (0.327")	n/a	10.0mm (0.394")	9.5mm (0.37")	9.5mm (0.37")	44.0mm (1.73")	47.6mm (1.89")	50 swabs per box 10 boxes per case



TX730 is a 100% nylon swab with micro bristles that grab particles or remove liquids by capillary forces. Ideal for tight spaces. Packaged in a silicone-free, amide-free bag.

Part No. TX730	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	0.8mm (0.031")	0.8mm (0.031")	15.5mm (0.610")	2.8mm (0.110")	2.8mm (0.110")	152.5mm (6.004")	168.0mm (6.614")	50 swabs per bag 10 bags per case



Sterile Available
All Texwipe swabs available in sterile upon request.

General Purpose Swabs

Texwipe's General Purpose Swab Series is designed for those customers needing a swab in industrial or general purpose applications. 100% polypropylene handle provides excellent chemical compatibility. Complete thermal bond construction eliminates adhesive contamination.

Packaged in a reclosable, silicone-free bag. Available in Polyester Honeycomb and Foam.



Reclosable, silicone-free bag

Polyester Honeycomb

Texwipe's General Purpose Polyester Swab Series, made with 100% polyester, was engineered with a honeycomb-type cell designed to maximize absorbency and capture of large particles and dirt in less critical areas.

Part No. TX801	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	13.5mm (0.531")	3.5mm (0.138")	25.2mm (0.992")	5.2mm (0.205")	3.0mm (0.118")	102.3mm (4.028")	127.5mm (5.020")	100 swabs per bag 6 bags per case



Part No. TX802	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.5mm (0.138")	3.0mm (0.118")	9.2mm (0.362")	2.5mm (0.098")	2.5mm (0.098")	60.0mm (2.362")	69.2mm (2.724")	500 swabs per bag 5 bags per case



Sterile Available
All Texwipe swabs available in sterile upon request.

Texwipe is certified to ISO 9001 and ISO 13485.



General Purpose Swabs



Reclosable, silicone-free bag

Foam

Texwipe's General Purpose Foam Swab Series, made with 100% polyurethane foam, was engineered for applying and absorbing adhesives, greases and liquids in less critical areas.


Part No. TX803	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.6mm (0.142")	3.5mm (0.138")	12.0mm (0.472")	3.0mm (0.118")	3.0mm (0.118")	59.0mm (2.323")	71.0mm (2.795")	500 swabs per bag 3 bags per case



Part No. TX804	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	3.5mm (0.138")	3.5mm (0.138")	15.5mm (0.610")	2.4mm (0.094")	2.4mm (0.094")	51.0mm (2.008")	66.5mm (2.619")	500 swabs per bag 3 bags per case



Part No. TX805	Head Width	Head Thickness	Head Length	Handle Width	Handle Thickness	Handle Length	Total Swab Length	Packaging
	19.7mm (0.776")	9.4mm (0.370")	26.4mm (1.039")	6.6mm (0.260")	3.2mm (0.126")	105.8mm (4.165")	132.2mm (5.205")	80 swabs per bag 5 bags per case

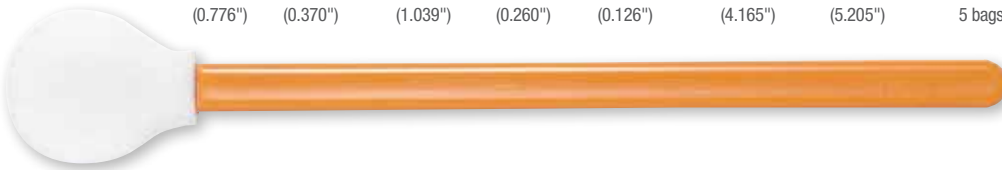

Sterile Available
All Texwipe swabs available in sterile upon request.



Swab Solvent Compatibility

Solvent	Head Substrate		Handle Material	
	Foam	Polyester	Polypropylene	Stat-Rite®
Acetic Acid	U	U		
Acetone	S	U	U	D
Acetonitrile	U	U	U	D
Ammonia, Anhydrous	S	U		
Aniline	S	U		
Benzene	S	U	U	D
Carbon Tetrachloride	S	U		
Chloroform	S	U		
Chlorophenol	D	D		
Chromic Acid	D	D		
Cyclohexane	S	U		
Dimethylformamide	D	S		
Dimethylsulfoxide	D	S		
Ethyl Acetate	S	U		
Ethyl Alcohol	U	U		
Ethyl Chloride	S	U		
Formic Acid	D	U		
Hexane	U	U	U	U
Hydrochloric Acid 48%	D	D	U	U
Hydrogen Peroxide	D	U		
Isopropyl Alcohol	U	U	U	U
Methanol	S	U	U	U
Methyl Ethyl Ketone	S	U		
Methylene Chloride	S	U	U	D
Nitric Acid 100%	D	D	U	D
Nitrobenzene	S	U		
Perchloroethylene	U	U		
Phosphoric Acid 85%	D	U		
Sodium Hydroxide 36%	S	U		
Sulfuric Acid 98%	D	D		
Toluene	S	U	U	D
Trichloroethylene	S	U		
Vinyl Acetate	S	U		
Water	U	U	U	U
Xylene	S	U	U	D

Compatibility is determined at room temperature.

Key: D – Degrades; S – Swells, U – Unaffected

Swab Index

STX705P	Pg. 2
STX705W	Pg. 2
STX708A	Pg. 2
STX712A	Pg. 2
STX763	Pg. 2
STX764	Pg. 2
STX705PT	Pg.3
STX764T	Pg. 3
TX700B	Pg.10
TX705	Pg.10
TX706A	Pg.7
TX707A	Pg.7
TX708A	Pg.7
TX709A	Pg.7
TX710A	Pg.7
TX712A	Pg. 7
TX714A	Pg. 6
TX714K	Pg. 4
TX714MD	Pg. 5
TX715	Pg. 4
TX716	Pg. 4
TX720B	Pg. 10
TX726	Pg. 10
TX730	Pg. 10
TX740B	Pg. 8
TX741B	Pg. 8
TX742B	Pg. 8
TX743B	Pg.6
TX750E	Pg. 9
TX751B	Pg. 8
TX752B	Pg. 8
TX753E	Pg. 9
TX754B	Pg. 6
TX757B	Pg. 8
TX757E	Pg. 9
TX758B	Pg. 6
TX758E	Pg. 9
TX759B	Pg. 6
TX761K	Pg. 4
TX758MD	Pg. 5
TX761	Pg. 6
TX761D	Pg. 9
TX761MD	Pg. 5
TX762	Pg. 6
TX769E	Pg. 9
TX801	Pg. 11
TX802	Pg. 11
TX803	Pg. 12
TX804	Pg. 12
TX805	Pg. 12

Sterile Available


All Texwipe swabs available in sterile upon request.

Texwipe is certified to ISO 9001 and ISO 13485.





Wipers

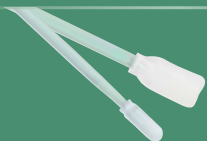

Swabs



Mops

Disinfectants



Solutions

Stationery



Mats

Sterile





Texwipe Americas +1 336 996 7046
Texwipe Europe +31 88 1307 410
Texwipe Asia +63 049 530 3390
Texwipe China +65 6468 9433

www.texwipe.com

Rev.01-09/15